

SENNEY

Collagen Light

Luxury Beauty Instrument

Quick Start Guide




SENNEY

COLLAGEN LIGHT LUXURY BEAUTY INSTRUMENT

Product Introduction

SENNEY COLLAGEN LIGHT LUXURY BEAUTY INSTRUMENT

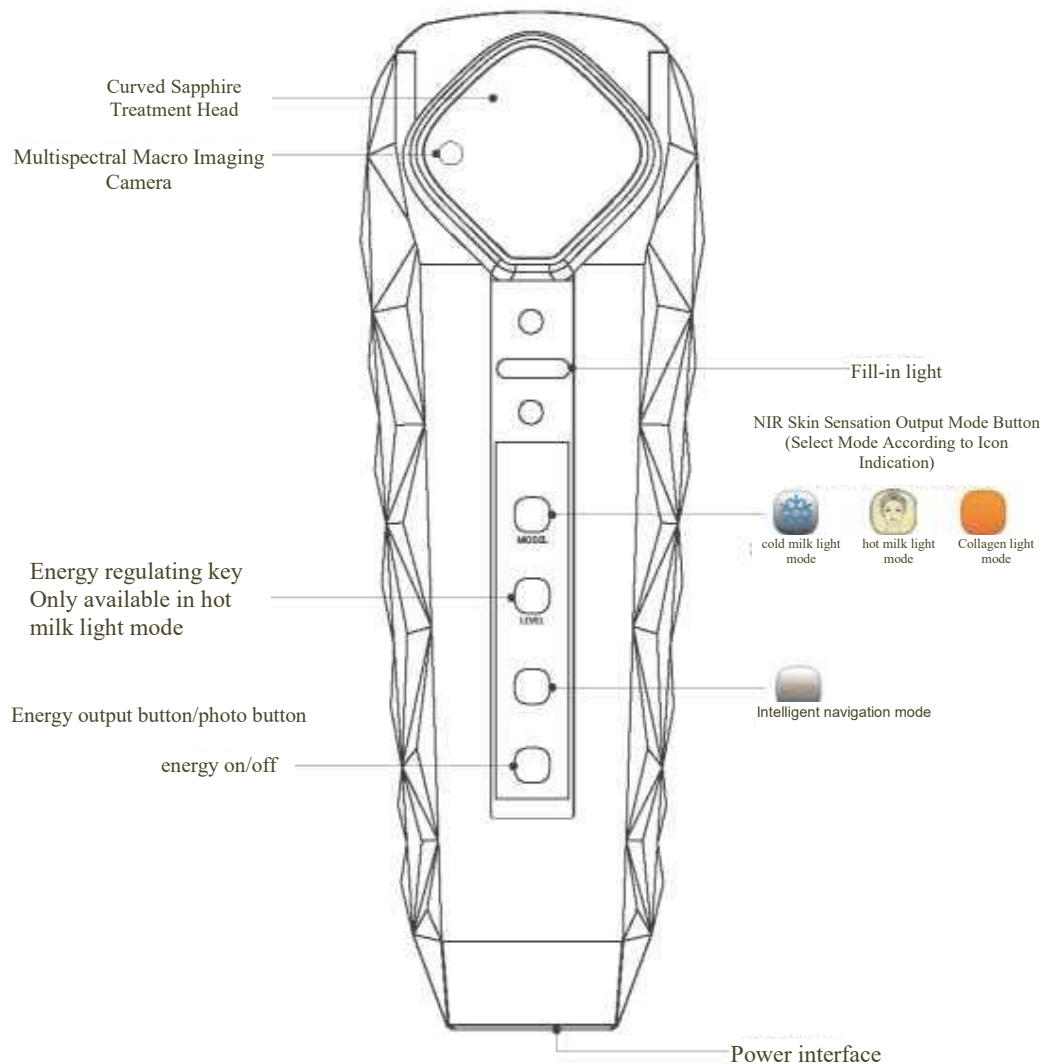

IMPORTANT NOTICE!!

To give you a better experience with SENNEY Milk Collagen beauty instrument, please be sure to read this manual carefully before your operation.

(The instrument display button adopts sliding touch operation and non-pressing touch)


Operating processes



❶ Clean your face with water after removing makeup. Be sure to remove contact lenses and accessories before use. Remove makeup and make your face clean.



❷ Long press the power button for 3 seconds to turn it on; the nursing head light flashes and the macro picture holding head self-light will be steady on. Then, the startup is completed.

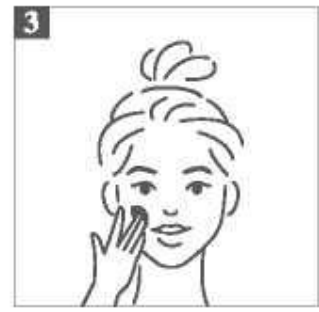


❸ Apply SENNEY Aloe Vera Moisturizing Gel for 1MM evenly to the operating area



❹ Choose corresponding nursing modes according to your own needs. There are two kinds of light for you to choose from.



❺ Hold the instrument and make its head close to your skin; draw it from the inside out and ensure its horizontal slide. Slide each part for more than 3 times.

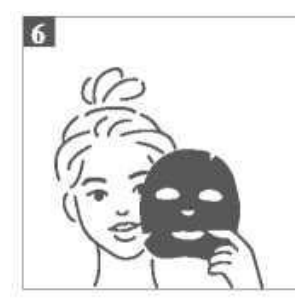

❻ After the instrument is used, turn off the power and wipe the excess gel off the face with wet wipes. Apply SENNEY Bosin moisturizing and firming mask. Wipe the head with wet wipes and finish the treatment. Pay attention to daily sun protection.



Cold Milk Light Mode

(Soothing and Calming)

Full-Screen Sapphire Cooling

Neutralize Burning Sensation, Soothe and Calm Sensitivity

5 Minutes per Use, Can Be Used Daily



Warm Milk Light Mode

(Brightens Skin Tone)

Adjust Energy Intensity Based on Skin's Heat Tolerance

Use This Mode Only After Becoming Familiar with the Device

15 Minutes per Use, Use Every 2 Days



Collagen Light Mode

(Stimulate Collagen Production)

Maximum Energy Output, Recommended for Mature Skin, Users with Good Skin Tolerance and Anti-Aging Needs

15 Minutes per Use, Twice a Week



Intelligent Navigation Mode

(Comprehensive Care)

Automatically Adjusts Energy Output Based on Skin Condition

One-Touch Comprehensive Skin Care

15 Minutes per Use, Twice a Week

Notice

Adaptation Test


Before using this product and the accompanying gel for the first time, please conduct a skin adaptation test.

Operations steps are as follows:

1. Clean the inner side of your arm with a facial cleanser and dry it with a towel.
2. Apply an appropriate amount of gel to the inner side of your arm, then use the Intelligent Navigation Mode for 3 minutes.

If the tested area shows no skin abnormalities (such as prolonged redness, burning, or overheating) within 24 hours, you may proceed to use the product. If any abnormalities such as rashes or itching occur within 24 hours before or after the test, discontinue the test and refrain from using the product.

Forbidden Area

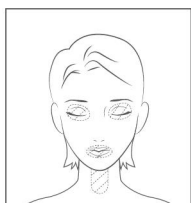

During facial care, please do not use this product to the following shaded areas:

Areas within the eye socket, including the eyeball and upper eyelid;

The throat, Adam's apple area, and lymph nodes

Warm Tip

Please apply the gel before activating the working mode; otherwise, sapphire on the treatment head may overheat.

Ensure that the device's contact head is fully in contact with the skin during use.

Do not stay on the same area for repeated use; always keep the device moving across the skin.

Do not rinse the device with water, and avoid wiping it with corrosive liquids such as alcohol, alkaline detergents, or nail polish remover.

Be sure to use the device with the dedicated gel provided with this product.

Safety instructions

Precautions and safety instructions!

Remember: This is a home milk collagen photocosmetic device that uses near-infrared light.

Please observe the following precautions:

The instrument can only be used for what is described in this user's manual.

Turn on the instrument only during operation.

Do not use the instrument if it is disassembled or broken.

Do not use this device near flammable liquids (e.g., alcohol, acetone).

Please read the product manual carefully before operation.

Please keep the instrument out of the reach of children.

Do not place, store, or use this device near water sources (such as sinks or bathtubs).

Usage taboo!

If you meet any of the following criteria, please do not use the SENNEY Milk Collagen Beauty Instrument.

Have an implantable pacemaker, suffer from arrhythmia, or any other known serious heart condition.

The operating area has implantable metal device or perforated with a decorative ring.

Any kind of suspicious lesion in the operating area.

During pregnancy or within 100 days of lactation.

Under 18 years old.

Have an impaired immune system caused by an immunosuppressive disease, such as AIDS, or are taking immunosuppressive drugs.

Have cancer or had cancer surgery within a year, or have had skin cancer.

Have a blood clotting disorder or diabetes.

Take blood thinning medications (anticoagulants).

Have any disease caused by heat, such as recurrent herpes simplex virus (within the operating range).

Any skin diseases such as sores, psoriasis, eczema, rashes, etc. within the scope of operation.

Any surgery, such as laser skin replacement, chemical peeling, or ongoing wound healing within 3 months of operation.

Plastic permanent fillers were implanted within the operating range.

Currently used medications, herbal preparations, vitamins or food supplements may cause brittle or dry skin.

If you suspect that you are suffering from any other diseases or medications you may be taking, which may cause adverse effects when you use SENNEY Milk Collagen Beauty Instrument, please consult a professional doctor.

Warning!

- Do not use SENNEY Milk Collagen Light Beauty Instrument on the upper eyelid or genital area.
- SENNEY Aloe Vera Moisturizing Gel is developed for the use of the instrument, so you will gain better effects when it is used together with the instrument.
- SENNEY Bossin Moisturizing and Firming Mask is developed for skin care after you use the instrument, and it is recommended to use them together.

Expected Improvements

Immediate Improvement Expectations

After the treatment, you will immediately feel a warm sensation in the treated area. Additionally, you will notice that your skin is firmer, smoother, and more radiant.

Long-Term Improvement Expectations

Your skin is expected to gradually improve within one month of use as new, healthy collagen continues to form. You will notice a significant reduction in visible wrinkles, a noticeable decrease in spots, and skin that is more elastic and radiant. After completing the full two-month usage period (using the product every other day), it is recommended to maintain the results with follow-up treatments 2-3 times per week. For consumers with sensitive skin or thin cuticle, it is advised to use the product no more than twice a week, starting with the Cold Milk Light or Warm Milk Light mode at a low intensity.

Possible side effects


Since the working mechanism of the SENNEYSENEY Milk Collagen Light Beauty Instrument is based on infrared light and thermal energy, the skin in the treated area may feel redness and a warm sensation after use. This simply indicates that the necessary thermal energy and infrared light waves have been applied to achieve optimal results. Such redness and warmth will typically subside shortly after the treatment. If the redness and warmth persist for more than 24 hours, please discontinue using the device.

! Please use the Milk Light high-gear mode and the Collagen Light mode. Due to significant individual differences in skin types, some users may experience discomfort when using high-intensity settings. If discomfort occurs, please switch to a lower intensity setting.

* The results of use depend on a variety of factors, and individual outcomes may vary. Therefore, the manufacturer and distributor do not provide any promises or guarantees regarding the effects of use.

Items included in the beauty kit

- **SENNEY product host**
- **Power adapter**
- **Product manual**
- **User Guide**

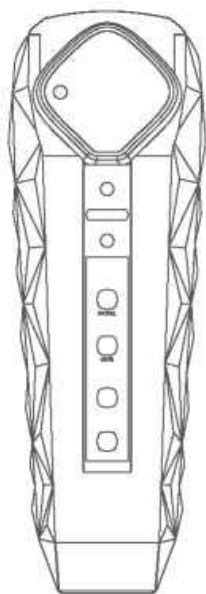

SENNEY product host



Power adapter



Product manual



User Guide

Basic product parameters

Product Name: NIR-PRO

Spectral Wavelength: 720-1900nm

Input Voltage: 24V Current: 1.5A

Dual Use: Rechargeable and Plug-in

Guarantee clause

1. Replace it only within a year

From the date of receipt of the order, the instrument with quality problem caused by non-human reasons will only be changed and not repaired.

2. The instrument can be returned within 7 days without reason if the secondary sale is not affected.

Within 7 days from the date of receipt of the order, if the product is not unsealed and does not affect the secondary sales, you can return or replace the product with no reason. If the product has been unsealed within 7 days after receipt, it is regarded as having been used, which will affect the secondary sales. If so, some depreciation fee will be collected for your return or replacement.

3. Related giveaways

If the main product is returned after sale, relevant giveaways would be returned together. Please ensure that the giveaways are complete and do not affect the secondary sales (if it has been unwrapped, 10% of the value of the gift will be charged).

Warranty certificate (any one of the following is acceptable):

1. Purchase Order Number
2. Purchase Invoice
3. Warranty Manual

QUALIFIED CERTIFICATE	
Product Name:	SENNEY Milk Collagen Light
Product Model:	Pro
Production Date:	
Compliance Standard:	GB4706.1-2005 GB/T36419-2018
Inspector:	
Manufacturer:	Shenzhen Pucci Technology Co., Ltd.
Factory Address:	503A, Building A2, Fuhai Industrial Zone, Fuyong Community, Fuyong Street, Bao'an District, Shenzhen

QC PASS inspection 01 (Sealed)



SENNEY