

SENNEY

SENNEY Precise Collagen Light—Instruction Manual



◀ Scan code to add SENNEY product experts



Whole course guidance



Exclusive benefit



2-year warranty



VIP

SENNEY

Thank you for purchasing SENNEY Precision Collagen Light

To ensure optimal performance and safety, please read this manual carefully before using the product and follow the instructions provided strictly.

After reading this manual, please keep this manual properly for future reference.

If the product is transferred to another user, please ensure that this manual is included.

Instructions for Use

Dear customer, thank you for purchasing the product SENNEY Precise Collagen Light of Shenzhen Pucci Technology Co., Ltd. Before using this product, please read and follow the contents of this instruction manual carefully so that you can use this product correctly and effectively. If you have any questions during the use process, you can consult the customer service on our company's e-commerce platform. We will wholeheartedly provide services to you.

Declaration

Our company shall only bear certain responsibilities for the safety, reliability and functionality of the equipment when all the following operations are complied with. Otherwise, our company shall not be liable.

The maintenance and adjustment of this machine shall be completed by the professional personnel of our company. If you disassemble the machine privately, after-sales service will not be provided.

Our product can comply with the dual standard certifications of national GB/T 36419-2018 and GB 4706.1-2005.

Please use the accessories provided with this product.

Product and Manufacturer Information

Product Information

Product Name: SENNEY Precise Collagen Light

Model: S2

DC Input Voltage: 24V

Current: 3A

Rated Power: 72W

Production date: Refer to Certificate of Conformity.



Manufacturer Information

Manufacturer: Shenzhen Pucci Technology Co., Ltd.

Factory Address: 503A, Building A2, Fuhai Industrial Zone, Fuyong Community, Fuyong Street, Bao'an District, Shenzhen

After-sales Service Provider: Shenzhen Pucci Technology Co., Ltd.

Address: 503A, Building A2, Fuhai Industrial Zone, Fuyong Community, Fuyong Street, Bao'an, Shenzhen

Adaptability Test

Before using the product and the supporting gel for the first time, please conduct a skin adaptability test.

Operation steps are as follows:

1. Clean the inner side of the arm with a facial cleanser product and dry the moisture with a towel.
2. Apply some gel to inner side of the arm and use the dual-effect collagen mode for 3 minutes.

If there is no skin abnormality (persistent redness that does not fade, burning, etc.) at the tested part within 12 hours, this product can be used. If abnormal conditions (rashes and itching) occur within 24 hours before and after the test, please stop the test and do not use this product.

Friendly Reminder

Please apply the gel to your face before starting the working mode.

Ensure the device head is in full contact with your skin during use.

Do not stay on the same area for more than 3 seconds; keep the device moving across the skin.

Do not rinse the device directly with water, and avoid using corrosive liquids such as alcohol, alkaline detergents, or nail polish remover to clean the device.

Always use the specialized gel provided with this product.

Precautions and Safety Instructions

Note: This is a home-use beauty device that utilizes near-infrared light, ultra-high-frequency electromagnetic waves, microcurrent, and 755NM precision semiconductor technology.

Please adhere to the following precautions:

The device is only intended for the purposes described in this user manual.

Turn on the device only during operation.

Do not use the device if it has been disassembled or damaged.

Before operation, carefully read the product instruction manual.

Keep the device out of reach of children.

Do not place, store, or use the device near water sources.

Usage Contraindications!

Please do not use SENNEY Precise Collagen Light if you meet any of the following conditions:

Having an implanted pacemaker, suffering from arrhythmia or other known severe heart diseases.

Having any implanted metal device or piercing with a ring in the operation area.

Having any form of suspicious lesion in the operation area.

During pregnancy or within 100 days of the lactation period.

Being under 18 years old.

Having an impaired immune system caused by immunosuppressive diseases, such as AIDS, or being on immunosuppressive medications.

Having cancer, having undergone cancer surgery within one year or having had skin cancer past.

Suffering from coagulation disorders or diabetes.

Taking blood-thinning medications (anticoagulants).

Having any disease triggered by thermal stimuli, such as recurrent herpes simplex virus (within the operation range). Having any skin diseases (sores, psoriasis, eczema, rashes, etc).

Having undergone any surgery, such as laser skin resurfacing, chemical peeling within 3 months in the operation range, or having wounds that are still healing, etc. Having implanted permanent cosmetic fillers in the operation range.

Using medications, herbal preparations, vitamins, or food supplements that may cause skin fragile or dry.

If you suspect that you have any other diseases or are taking medications that may have an adverse effect, please consult a professional doctor.

Warning!

Do not use the SENNEY Precise Collagen Light on the upper eyelids or the genital area.

The SENNEY Aloe Vera Moisturizing Gel is developed for the optimal performance of the device, and using it in combination will achieve better results.

The SENNEY Pro-Xylane Luxurious Moisturizing and Firming Mask is developed for skin care after using the device, and it is recommended to use it in combination.

Product Diagram

Precise Collagen Light


Key function

- 1. Press and hold the power button to turn the device on or off.
- 2. Press the power button twice in quick succession to switch the mode.
- 3. Press the power button once to adjust the gear level.

Key function



Attention

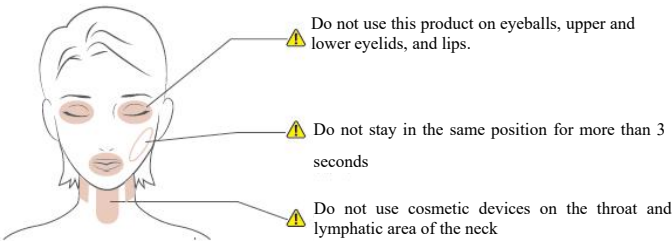

Product composition




1

SENNEY precision collagen light beauty tool x1



2

Operation manual x1

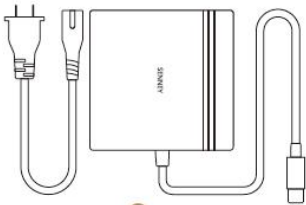

3

Power adapter x1

Description of functions

Product name		Precise Collagen Light S2
Function collocation (Fixed 15 minutes per mode)	V-FACE V-face Mode	Ultra-high frequency electromagnetic wave+ice compress + microcurrent
	F-COLLAGEN Dual-effect Collagen Mode A	NIR collagen light+ultra-high frequency electromagnetic wave +ice compress
	NIR Skin Rejuvenation and Brightening Mode	NIR collagen light + ice compress + microcurrent
	R-CEORO Ice Sensation Repair Mode	NIR collagen light + ice compress
Instrument material	Light-emitting Mirror Surface	Light-transmitting crystal
	Electrode Head	Titanium alloy/Medical grade 316 stainless steel
	Shell	Environmentally friendly engineering plastic/Aluminum alloy

V-FACE V-Lift Mode
High-Frequency Electromagnetic Waves
+ Microcurrent + Ice Therapy
Sliding Lifting Technique



Paired with **SENNEY Aloe Hydrating Gel**

F-COLLAGEN Dual Collagen Mode
NIR Collagen Light + High-Frequency
Electromagnetic Waves + Ice Therapy
Sliding Lifting Technique



Paired with **SENNEY Aloe Moisturizing Hydra-Gel**

NIR Skin Refining and Brightening Mode
NIR Collagen Light + Microcurrent + Ice Therapy
Sliding Lifting Technique



AA Paired with **SENNEY Aloe Moisturizing Hydra-Gel**.

R-CEORO Ice Cooling Repair Mode
NIR Collagen Light + Ice Therapy
Sliding Lifting Technique



Usage Recommendations for the Precise Light Mode

Due to individual differences, when using the precise light mode, you need to select an appropriate usage method according to your own skin condition. Do not overuse it in pursuit of quick results.

Recommended Usage Method (for Beginners)

Wear goggles and slide the device closely against the skin. Slide it 3 times on each area.

Advanced Usage Method (after using the recommended method for more than 2 weeks)

Irradiate the nursing area with the precise light, and stay for 1-2 seconds. When you feel a hot and burning sensation on the skin, you can change to the next area. Repeat the nursing process 3 times for each area.

After the nursing is completed, use the ice compress and repair mode of the Precise Collagen Light to ice compress the nursing area.

Possible Skin Reactions

In areas with severe spots and pigmentation, scabbing, peeling, or blistering may occur. You can apply antibacterial drugs such as growth factors and erythromycin eye ointment, and suspend the use of the device. After about 3 days when the skin stabilizes, you can continue to use it according to the recommended method.

Please note the following points before using the precise light mode:

1. If you have undergone picosecond laser treatment recently, please wait for 30 days before using this device.
2. Do not use it if your skin is unstable or has wounds.
3. People under 16 years old, pregnant women, and those with light-sensitive skin should not use it.
4. It is prohibited for those with severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin diseases, burned skin, and anticoagulation disorders.
5. For the areas where botulinum toxin or hyaluronic acid has been injected, please wait for 30 days before using. Please consult a professional regarding whether other subcutaneous fillers conflict with this product.
6. Avoid using it on sarcomas, moles, and nevi during use.
7. Do not use it to irradiate the eye area or directly aim the device at the eyes. It can only be used around the eyes where you can feel the subcutaneous bones. Otherwise, it may cause serious eye damage, such as permanent or temporary vision loss or blindness.
8. People under 16 years old, pregnant women, and those with light-sensitive skin should not use it. It is prohibited for those with severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin diseases, burned skin, and anticoagulation disorders.
9. Do not use it if you have abnormal collagen metabolism, including a history of scars and poor wound healing.
10. Do not use it if you have vascular diseases, such as varicose veins in the irradiated area.
11. It is prohibited to use it on inflamed areas, mucous membranes, and other sensitive parts.

This list is not exhaustive. If you are unsure whether you can use this product, it is recommended that you consult your doctor.

Nursing Recommendations in Normal Skin Condition

The normal skin condition refers to a state where there are no major skin problems such as allergies, acne breakouts, dryness, or peeling, and the skin is clean and has a uniform color.

Establish Tolerance within 1-2 Weeks

Use it at a low gear level for 2 weeks to help the skin gradually adapt to the temperature and the sensation of electric current brought by the photoelectric energy. After the skin gets adapted, you can gradually increase the gear intensity.

Advanced Anti-aging after 2 Weeks

After the skin has adapted to the photoelectric energy of the device, you can try using a higher gear intensity to achieve better results.

Precise Collagen Light S2

	Monday	Tuesday	Wednesday	Thursday	Friday	Saturday	Sunday
V-FACE- mode	√		√		√		√
F-COLLAGEN mode	√		√		√		√
NIR mode	√		√		√		√
R-CEPRO mode	√	√	√	√	√	√	√

Warranty and Warranty Certificate

Mainland China Region Our company provides 1-year quality warranty services for customers purchasing this product in Mainland China. The warranty period is calculated from the goods receipt date. For electronic transactions, it shall be based on logistics receipt date; for physical transactions, it shall be based on physical delivery date.

Our company complies with relevant laws and regulations, including *E-commerce Law of the People's Republic of China*, the *Product Quality Law of the People's Republic of China*, and the *Law of the People's Republic of China on the Protection of Consumers' Rights and Interests*. For offline physical transactions, if a quality problem occurs with the purchased product within 7 days (inclusive) from the purchase date, customers can contact the seller for return or exchange with an invoice or a similar valid certificate; if a quality problem occurs within 15 days (including) from the purchase day, customers can contact the seller for exchange with an invoice or a similar valid certificate. For online electronic transactions, customers can return the goods without no reason in 7 days from goods are received. The returned goods should be in good condition and not affect secondary sales. The freight for 7-day no-reason return shall be paid by customers, except for those caused by product quality. If you purchase this product outside Mainland China, the product's warranty period, return and exchange policy, maintenance and other related matters will comply with laws and regulations of that country or region. For specific matters, please contact the authorized dealer in your location for assistance.

Please request a purchase certificate (an invoice) when purchasing this product and keep it properly. During the warranty period, if the product's warranty period cannot be determined due to the inability to present an invoice or a valid purchase certificate, the production date on the product's body production label will be used as the starting time of the warranty period.

Please note that the following situations will not be eligible for free warranty service during the warranty period:

- Artificial damage caused by improper use, such as connecting to an inappropriate power source, using inappropriate accessories, not using it according to the instruction manual, incorrect use, or damage caused by negligence, etc.
- Maintenance or modification not authorized or approved by our company.
- Damage caused by other force majeure factors (such as natural disasters, abnormal voltage, etc.).
- Product aging and wear caused by normal use, provided that it does not affect the normal use of the product.

The free service does not include consumable spare parts and accessories that need to be replaced frequently, such as the power adapter, etc.; the free service does not include transportation fees, etc.

SENNEY

After-sales Three Guarantees Service

- Within 7 days starting from the next day after the signing for receipt, if a performance failure occurs and is confirmed by the after-sales inspection of SENNEY, you can enjoy the free return and exchange service;
- Within 3 months starting from the next day after the signing for receipt, if a performance failure occurs and is confirmed by the after-sales inspection of SENNEY, you can enjoy the free replacement and repair service;
- Within 12 months starting from the next day after the signing for receipt, if a performance failure occurs and is confirmed by the after-sales inspection of SENNEY, you can enjoy the free repair service.

Non-Worry-Free After-sales Regulations

1. Things have exceeded the validity period of the three guarantees/do not conform to *Performance Failure Standard*.
2. Unauthorized maintenance, tearing or altering the label or anti-counterfeiting mark.
3. Human factors cause performance failures of this product and its accessories.

Performance Failure Standard

1. Abnormal power on and off.
2. Abnormal charging.
3. Unable to work normally.

After-sales Warranty Certificate (Any one of the following can be used as a warranty certificate)

1. Purchase order number.
2. Purchase invoice.
3. Warranty manual.

SENNEY